

Tiga FX

RGK



The most important quality of the Tiga FX is its Q-lock system. Innovative, smart and revolutionary, you can use these to reinvent your rigid made to measure wheelchair to a collapsible, petite travel companion. Remove the cushion, flatten those sideguards, fold the front frame underneath the seat rail, the backrest takes the place of the cushion, just the quick release axle pins to go... get on with your journey.

Unique 'Q Lock' Folding System.
Image 1.

Q Lock folding frame and backrest.
Image 2.

"Fold flat" folding sideguards. **Image 3.**

Airtech upholstery system. **Image 4.**



Image 1



Image 2



Image 3

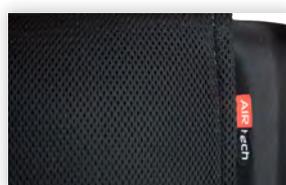


Image 4

Keeping you in mind. More options.

Quality products with the required accessories

SPECIFICATIONS

Frame Material	7020 Aluminum
Frame style	Rigid
Seat Width	10" - 18" / 25 cm - 46 cm
Seat Depth	10" - 18" / 25 cm - 46 cm
User Weight Capacity	125 kg / 275.6 lbs
Product Weight	From 8.5 kg / 18.7 lbs
Folding system	Q-lock folding system
Available colours	25 options



Rear wheels
Quick release

BENEFITS

For active people mainly motor injuries.

Discover our products



DESIGNED FOR OUR TERRITORIES

Our products are designed to meet the needs of each of our territories, offering quality and affordable prices.



FUNCTIONAL VERSALITY

We have a variety of products ranging from a wide range of manual and power wheelchairs and mobility scooter.

**For more information on obtaining this product
Contact your local Loh Medical representative**

Find your local representative at www.lohmedical.com   

ver_2025.11.11E

